



AVENIO NGS Oncology Assays

Comprehensive genomic profiling (CGP), targeted profiling, and surveillance & monitoring of solid tumors for clinical research.

For Research Use Only. Not for use in diagnostic procedures.



FOUNDATION
MEDICINE®




A versatile solution to meet your research needs

AVENIO NGS Oncology Assay kits provide researchers a versatile solution for tumor tissue and ctDNA analysis with an integrated end-to-end workflow. With accurate and reliable analysis of ultra low ctDNA detection as well as formalin-fixed paraffin-embedded (FFPE) tissue samples, AVENIO NGS Oncology Assays offer laboratories in-house tumor profiling, comprehensive genomic profiling (CGP), and surveillance & monitoring. Each kit includes reagents for DNA extraction, sample input QC (for tissue), library preparation, target enrichment, robust bioinformatics, secondary analysis and reporting (variant calls).

Select the assay that best meets your research needs

A portfolio of high-performance kit-based assays for in-house tumor profiling, comprehensive genomic profiling (CGP), and surveillance & monitoring for clinical research.

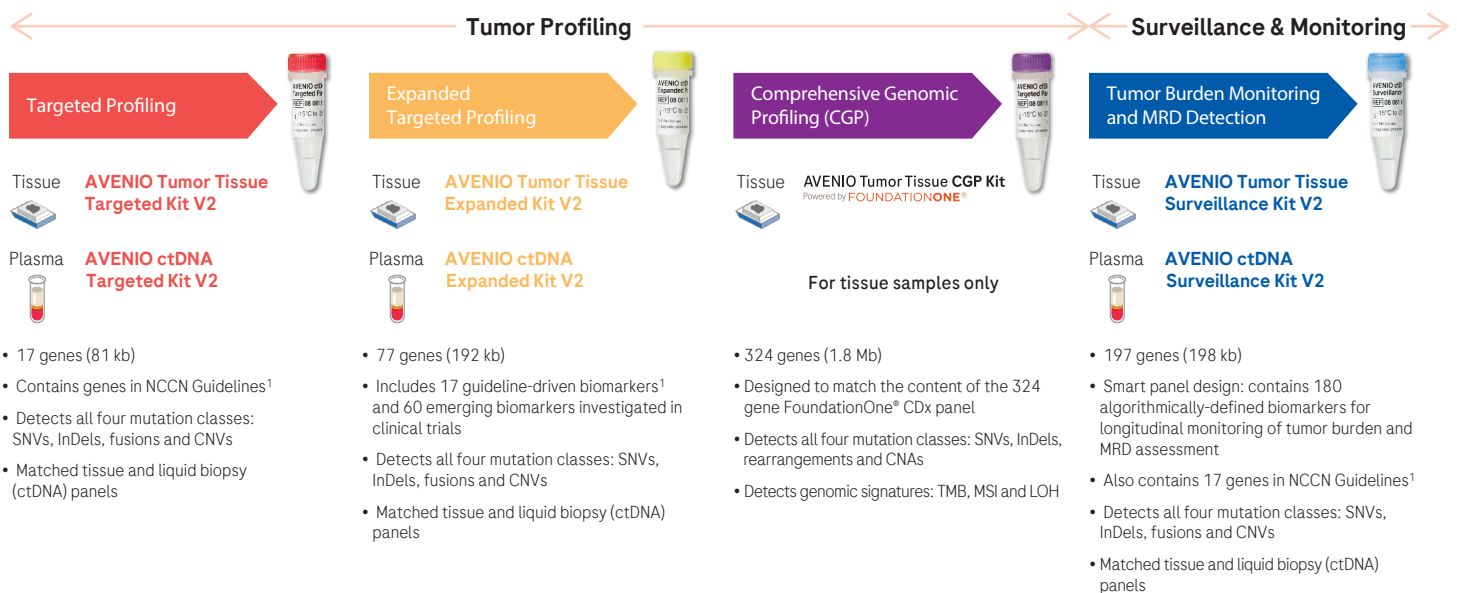

For Research Use Only. Not for use in diagnostic procedures.


¹ National Comprehensive Cancer Network. <http://www.nccn.org>. Accessed Aug 29, 2023.

Example research use cases

Example A

AVENIO Tumor Tissue CGP Kit
Powered by **FOUNDATIONONE**


Assess 324 cancer-relevant genes plus TMB, MSI and LOH in a single panel that is designed to match the genomic content of FoundationOne® CDx panel

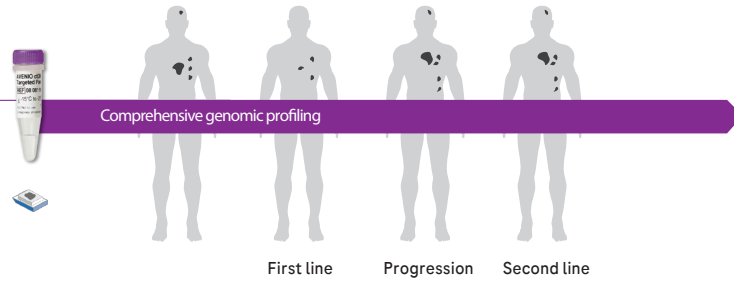

Example B

Targeted Kit


Identify biomarkers in NCCN Guidelines¹

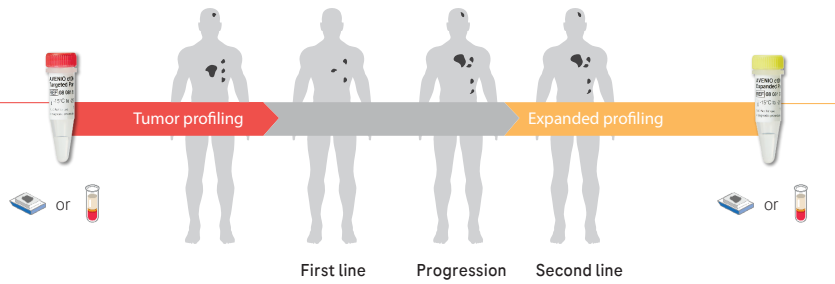

Example C

Surveillance Kit

Monitor tumor burden longitudinally in lung cancer and colorectal cancer

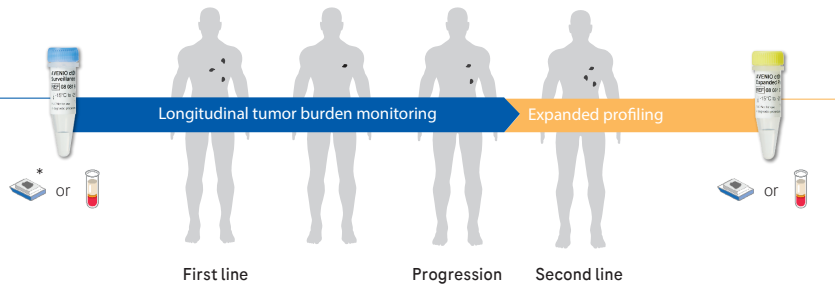


Expanded Kit

Identify biomarkers in NCCN Guidelines¹ and emerging biomarkers investigated in clinical trials

Example D

Surveillance Kit

Surveil for minimal residual disease in lung cancer and colorectal cancer

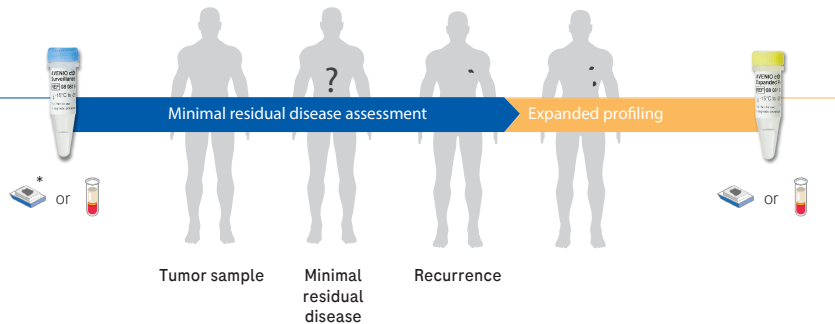

For Research Use Only. Not for use in diagnostic procedures.

* The ctDNA Surveillance Kit V2 may also be used if a tumor tissue sample is available.

¹ National Comprehensive Cancer Network. <http://www.nccn.org>. Accessed 29 August, 2023.

Performance you expect from Roche and Foundation Medicine


AVENIO Tumor Tissue CGP Kit: Analytical variant detection performance across genomic alterations and signatures

Libraries were prepared from 314 FFPE-derived DNA samples by the AVENIO Tumor Tissue CGP kit. For each alteration classification, the percentage of expected variants that were detected by the AVENIO Tumor Tissue CGP kit are shown. For genetic signatures, the percentage of expected samples detected as MSI, TMB, and LOH high assessment are shown.

Classification	Detected Variants/Signatures
Short Variants	98.2%
Rearrangements	90.5%
CNA	94.8%
MSI high	100%
TMB high	100%
LOH high	96.8%

Exceptional Performance as demonstrated by Key Sequencing Metrics

Libraries were prepared from 314 FFPE-derived DNA samples by the AVENIO Tumor Tissue CGP kit. Eight samples were sequenced per NextSeq 500 High-output flowcell. The graphs show sequencing QC metrics from the FoundationOne® Analysis Platform.

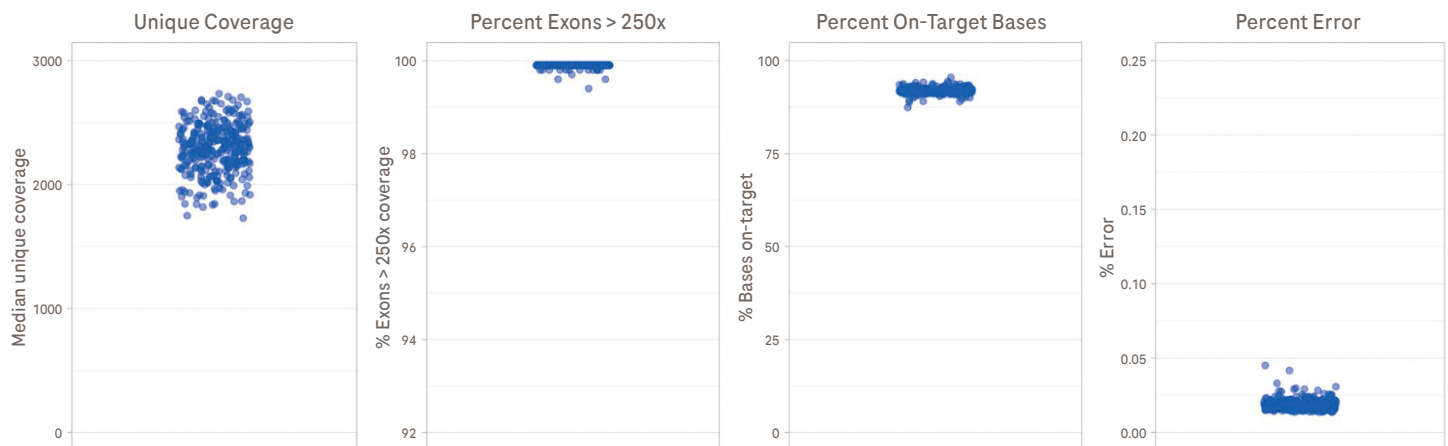

For Research Use Only. Not for use in diagnostic procedures.

Data on file with Roche.

FFPE, Formalin-fixed paraffin-embedded. QC, Quality control. CNA, copy number alterations. MSI, Microsatellite instability. TMB, Tumor mutational burden. LOH, Loss of heterozygosity.

An integrated solution for an end-to-end hybrid-capture workflow

The efficient, high-quality AVENIO workflow includes all sample prep reagents, input QC (tissue only), robust bioinformatics and secondary analysis all from single trusted source. It has been optimized to minimize hands on time, and deliver high-quality results in just 5 days, making it easy for you to obtain reliable genomic insights about solid tumors in your lab.


5 day workflow from DNA isolation to data analysis

Day 1	Day 2	Day 3	Day 4	Day 5
DNA Isolation	Library Generation	Target Enrichment	Sequencing	Analysis & Report

AVENIO Tumor Tissue Analysis Kits V2



AVENIO Tumor Tissue CGP Kit

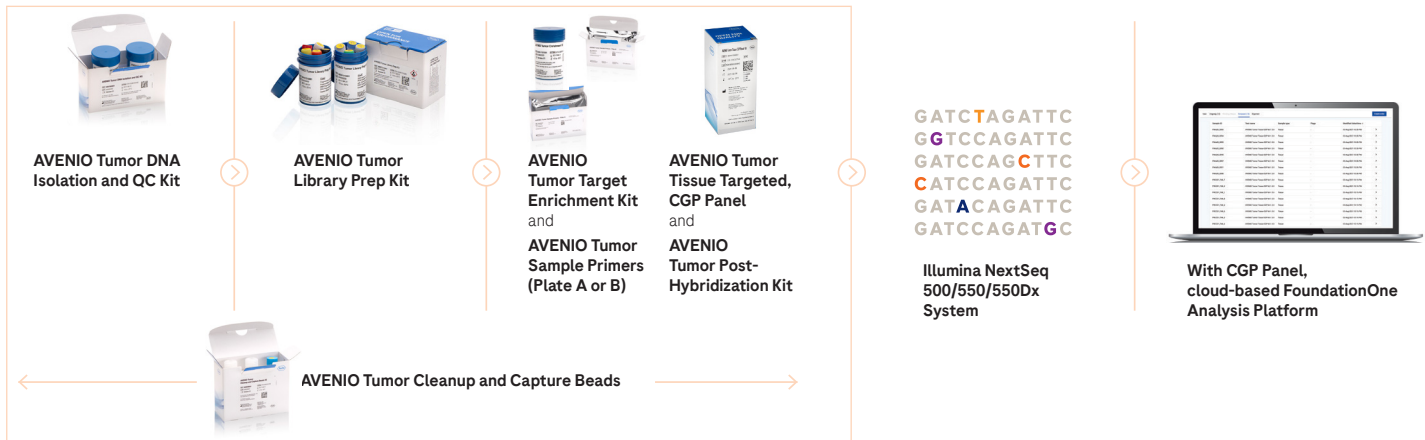

AVENIO ctDNA Analysis Kits V2



Exceptional performance¹

Rigorously optimized using thousands of samples, AVENIO assays deliver the level of performance you expect from Roche, a world leader and trusted partner in oncology.

AVENIO Tumor Tissue Analysis Kits V2

Mutant Allele Frequency/ Copy Number	SNVs		Indels		Fusions		CNVs	
	5%		5%		5%			
AVENIO Tumor Tissue Kit V2	Sensitivity	PPV	Sensitivity	PPV	Sensitivity	PPV	Sensitivity	PPV
Targeted	> 99%	> 99%	> 99%	> 99%	> 99%	> 99%	Sample dependent	> 99%
Expanded	> 99%	> 99%	> 99%	> 99%	> 99%	> 99%	Sample dependent	> 99%
Surveillance	> 99%	> 99%	> 99%	> 99%	> 99%	> 99%	Sample dependent	> 99%

Samples: FFPE tissue curls/sections, DNA input: ≥20 ng of FFPE DNA, total DNA amount for each sample determined by input QC.

Sensitivity and Positive Predictive Value (PPV) metrics based on observed product performance. Sensitivity and PPV performance reported per variant. Sensitivity was determined using commercially available reference samples containing verified mutations at the stated allele frequencies. SNV performance is panel wide. Indel, Fusion, and CNV performance based on whitelist variants. AVENIO Tumor Tissue Analysis Kits V2 achieve >99.99% per base specificity across each of the panels. Stated performance requires at least 20 million reads per sample for Targeted, Expanded, and Surveillance Kits V2. Sequencing performed on an Illumina NextSeq 500/550/550Dx instrument.

AVENIO ctDNA Analysis Kits V2

Mutant Allele Frequency/ Copy Number	SNVs		Indels		Fusions		CNVs**	
	0.5%*		1.0%*		1.0%			
AVENIO Tumor Tissue Kit V2	Sensitivity	PPV	Sensitivity	PPV	Sensitivity	PPV	Sensitivity	PPV
Targeted	> 99%	> 99%	> 99%	> 99%	> 99%	> 99%	Sample dependent	> 99%
Expanded	> 99%	> 95%	> 99%	> 99%	> 99%	> 99%	Sample dependent	> 99%
Surveillance	> 99%	> 99%	> 99%	> 99%	> 99%	> 99%	Sample dependent	> 99%

* Detects variants down to 0.1%

**Can detect ERBB2, EGFR and MET CNVs. Performance dependent on input sample type and input concentration (10 ng - 50 ng per reaction).

Performance samples: Seraseq ctDNA Complete Mutation Mixes (SeraCare), healthy donor cfDNA.

Sensitivity and Positive Predictive Value (PPV) metrics based on observed product performance. Sensitivity and PPV performance reported per variant. Sensitivity was determined using commercially available reference samples containing verified mutations at the stated allele frequencies. SNV performance data based on hotspot calls. The AVENIO ctDNA Analysis Kit V2 can achieve >99.99% per base specificity across each of the panels. Stated performance requires at least 40 million reads per sample for all AVENIO ctDNA Panels V2. Sequencing performed on an Illumina NextSeq 500/550/550Dx instrument.

For Research Use Only. Not for use in diagnostic procedures.

¹ Data on file with Roche.



OPEN FOR RESULTS

AVENIO Tumor DNA Isolation and QC Kit

REF: 08357293001	GTIN: 07613336154182
LOT: 00000000	2017-08-31
Version 01	-15 to -25°C

For Research Use Only. Not for use in diagnostic procedures.

For instructions and safety information, visit: <http://hangenp.org/roche.com/support>

Manufactured in Germany
Roche Molecular Systems, Inc.
480 Giltspur Street
Juelich, Germany
AVENIO is a trademark of Roche. © 2017 Roche Molecular Systems, Inc.

Distributed in USA
Roche Diagnostics
Indianapolis, IN 46202 USA
© 2017 Roche Diagnostics
48121 7819

Roche Diagnostics GmbH
68229 Mannheim, Germany
48121 7819

Bring the power of Roche and Foundation Medicine* into your lab

Contact your Roche representative today



Proven
Technology



Comprehensive
Portfolio



Trusted
Expertise

For Research Use Only. Not for use in diagnostic procedures.

*for the AVENIO Tumor Tissue CGP Kit only.

© 2023 Roche Sequencing Solutions, Inc. All rights reserved.

AVENIO is a registered trademark of Roche. FOUNDATIONONE is a registered trademark of Foundation Medicine, Inc. All other product names and trademarks are the property of their respective owners.



Bring our expertise in-house at sequencing.roche.com/aveniocgpkit or contact your local Roche representative for more information.

MC-09193